


70.00 mm

Dear customer, thanks for purchasing this product. Please read this user manual carefully and follow the correct instructions.



PROFESSIONAL DARK SPOT REMOVER

90.00 mm

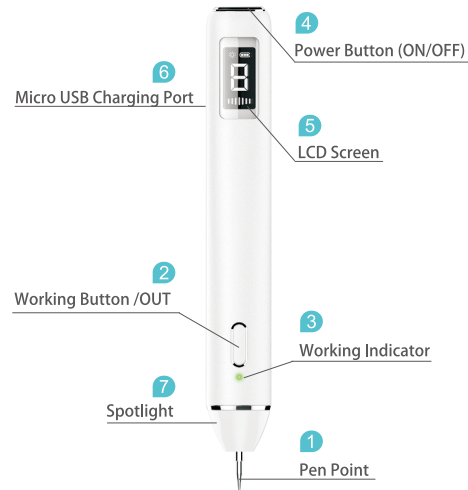

The Professional Dark Spot Remover adopts latest high technology, helps you to remove skin tags, dark spots, freckle, mole, wart etc effectively and instantly. Easy and Safe to use. With this beauty device, you will improve your skin and change your life!

Package Contents


- 1 * Skin Tag Remover
- 10* Fine Needles
- 1 * USB Cable
- 1 * User Manual
- 3 * Coarse Needles (including an extra needle guard)

Specifications


Material	ABS+Steel Needle	Input Current	5V/1A
Battery Capacity	800mAh 3.7V	Output Power	MAX5W
Product Size	182* 26* 26mm	Product Weight	75g
Package Size	200* 96* 36mm	Package Weight	200g

Features

- The 6th generation of German technology manufacture
- Microcomputer controlling and overdischarge protection, safe to use
- 10 adjustable levels meet different needs
- LCD screen shows the level and battery gauge
- Spotlight makes it convenient when work in dark places
- 10 replaceable fine needles and 3 coarse needles, sturdy and durable
- The battery icon will flash when the battery is low

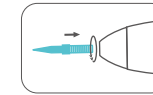

- 1 Thread coarse needle into the pen tip directly.

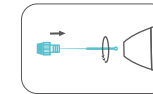

- 2 Insert fine needle into the needle guard first, then thread the needle guard into the pen tip.

Application Methods

1. Assemble the needle : disinfect the needle and skin with alcohol. Screw the coarse needle into the tip hole clockwise. Or insert fine needle into the needle guard and screw the needle guard into the tip hole clockwise.
2. Turn on/off : press the power button for 2s to turn on/off. The LCD screen shows the battery gauge and the level. If there is no operation for 10 mins, it will power off automatically.
3. Control gears: press the power button for one second to switch the gear from 1-10 to choose a decent level. Holding the pen tight while you press the OUT button which will turn on the green light. Gently scratch the pore or skin back and forth. The needle cannot stay in one place for too long and it will turn to sleep mode automatically if it keeps working more than two minutes.
4. LED light: while the pen is powered on, a quick double click of the power button will turn on the spotlight, and you can turn it off by double click it again.

Charging Methods

1. Use a 5V/1A adapter (prepared by customers) to connect the USB cable and the pen. The green indicator flashes.
2. The green light represents fully charged (about 2 hours).

Please do not use the device during charging. Besides, overcharge or charging when the battery is exhausted will shorten the life span of the battery.

Level Application (for reference only)

Level	Needle Type	Apply to
1-2	Fine	Fat granule
3-4	Fine	Small moles
	Coarse	Small spots
5-6	Fine	Small & Medium-sized moles
	Coarse	Small & Medium-sized spots
7-8	Fine	Mung bean-sized moles
	Coarse	Mung bean-sized spots and plaques
9-10	Fine	Large moles
	Coarse	Large spots, plaques, etc.

Caution: Users must start using from level 1, and select other levels according to the using results.

Notes

1. Before use the device on your skin, please practice on animal skin until you can operate it expertly.
2. You should consider applying partial anesthetic before the treatment according to your pain tolerance and your treating part. Feel comfortable for the smell of burning during the treatment.
3. Please do not scrub or clean the treating area with strong alkaline liquid. A simple care will be OK. Please do not eat spicy food and keep a healthy lifestyle during the healing period (1-3 months).
4. Please do not treat a large area (larger than 2-3cm) in one course of treatment as large wound increases the risk of infection and bleeding. It is recommended to divide into four or five courses of treatment.
5. When the fresh meat grows, it may feel itching. Please do not scratch the scab and it will peel off itself in about 20 to 30 days.
6. It's better to wear gloves during operating.
7. Please take the needle out timely after use and keep the needle in dry place after cleaning and disinfecting.
8. Please do not immerse this device into water for washing.
9. Please do not use the device while charging. It is recommended to charge every three months to protect the life span of battery when the device remains no-operation for a long time.
10. Please keep the device away from children, pets and irrelevant personnel.
11. The use of elders, children, pregnant women, people having heart diseases and people susceptible to scare must be under great caution. Please consult professionals if there is any doubts.

Practical Reference

Removing Dark Spots: you can see clearly that the surface of melanin turns charred after making the acusector contact the spot slightly and moving it back and forth. Then wipe the spot off with cotton swap. If the spot is light and tiny, treating for just one time will be enough. You can move back and forth for different times depending on the sizes of spots. Attention to the operation process, the experienter should keep still and the operator should hold the pen tightly without hand trembling to avoid burning and leaving any scar.

Removing Moles: it is very important to tell the size and depth of the moles when removing them. Some moles grow in the surface of skin, but some deep in the subcutaneous tissue. It is recommended to remove moles in small squares rather than prick the whole surface of the mole by one time. Besides avoiding inflammation, you'd better, according to doctor's advice, apply some drugs or nutrient cream conducive to cell growth to avoid skin irritation 2-3 days before the wound healing.

Removing Granulations: slightly contact the root of granulation with the acusector to burn the root of the derived granulation with laser. And the granulation will fall off automatically in a week or so. The removal of deep growth granulation can be divided into several times. Besides avoiding inflammation, you'd better, according to doctor's advice, apply some drugs or nutrient cream conducive to cell growth to avoid skin irritation 2-3 days before the wound healing.

Removing Sarcoids: slightly contact the root of sarcoid with the acusector to burn the root of the derived sarcoid with laser. And the sarcoid will fall off automatically in a week or so. The treatment of deep growth sarcoid can be divided into several times. Besides avoiding inflammation, you'd better, according to doctor's advice, apply some drugs or nutrient cream conducive to cell growth to avoid skin irritation 2-3 days before the wound healing.

Removing Tattoos: it is better to use coarse needles as the areas of removing tattoos are generally larger. Except costing more time, the method of removing tattoos is almost the same as that of removing dark spots. The removal of larger tattoos (more than 3cm) should be divided into several times. Please do not finish the removing in one time or it will cause bacterial infection and disastrous consequences. Users should be quite cautious. If severe inflammation happens, go to hospital as soon as possible.

90*70mm